



Pipettes

Our Leading Range of Products

Pipettes designed for micromanipulation, denuding and handling

OOCYTE RETRIEVAL | ANDROLOGY & IUI | FERTILIZATION | CULTURE | PGD & PGS | CRYOPRESERVATION | EMBRYO TRANSFER



CooperSurgical Companies

The Most Comprehensive Range of Micropipettes Available Today

Gives You Choices

Pipettes to suit any preference you may have for microtools for ICSI, biopsy, denudation and handling.

The Quality You Expect

Strict quality demands and rigorous QC ensures the highest possible quality and consistency of all products.

The Three Company Stories



In 1992 the first commercial high quality micropipets for IVF professionals was marketed by HUMAGEN. Product development of the micropipets has always been based on customers' needs and requirements. HUMAGEN micropipets are now marketed as ORIGIO micropipets – the quality and consistency of the products remain at the same high level.

Realizing the need for an easy tool for oocyte denudation The Stripper® was developed in the late 1990s. This revolutionizing product removed the need to prepare and pull your own denudation pipettes and was quickly an enormous success with ART professionals around the world.



Research Instruments was founded in 1962. As a global leader in micromanipulation equipment technology, the company moved into the fertility industry in the early 1990s. Working closely with pioneers in the field RI refined its micromanipulators to meet the needs of those perfecting ground breaking techniques such as ICSI.

Based on customer requests and as a natural extension of the micromanipulator and laser products, RI later marketed high quality micropipettes as well as highly user friendly denuding and handling pipettes.



The Pipette Company was founded in 1999 by a group of Australian scientists working in reproductive medicine. Strict principles of designing micropipettes for optimal clinical outcomes made the high quality pipettes a success from day one.

In 2000 the first micropipettes were shipped to embryologists and the popularity of the very carefully designed ICSI pipettes grew quickly.



CooperSurgical

CooperSurgical acquired ORIGIO in 2012, Research Instruments in late 2015 and The Pipette Company in early 2016.

Now, as CooperSurgical Companies we have an even stronger joint offering of high quality products that empower clinics to succeed.

Dedication

To our customers, to our products

Expertise

Knowledge of ART processes and products and decades of experience

Quality

Highest possible quality products that meet all standards for medical devices



CooperSurgical Companies



ICSI Pipettes

Choose from the largest range of standard ICSI pipettes ever offered - all carefully designed to meet your preferences and help you optimize outcomes.

Customer Favorites:

ORIGIO SlimLine Pipet - MIC-SLM-30

- Short taper for added stability when injecting
- Smaller inner diameter to minimize the impact on the oocyte



RI ICSI-Plus™ Spiked Injection Pipette - T5SI30

- Long taper for a more flexible tip
- Easy to use packaging - the TwistaPak™ opens like a lipstick making it simple to extract the pipette



TPC Sharp-Blunt Spike™ Pipette – LISR

- Extra long parallel taper, specifically designed for very fine fluid control and suction during use
- Sharp-Blunt Spike for easy penetration of the zona pellucida and gentle progression through the oolemma to minimize risk of tearing



Our variety of ICSI pipette designs

ORIGIO Pipet

RI Pipette

TPC Pipette

Holding Pipettes


A variety of carefully designed holding pipettes are available that ensure stability of your oocyte during ICSI.

Large and extra large holding pipettes are designed for extra stability of embryos during biopsy and other micromanipulation procedures.



Biopsy Pipettes

Our range of pipettes designed for biopsies are available in different sizes and design for polar body, blastomere or trophectoderm biopsy procedures.

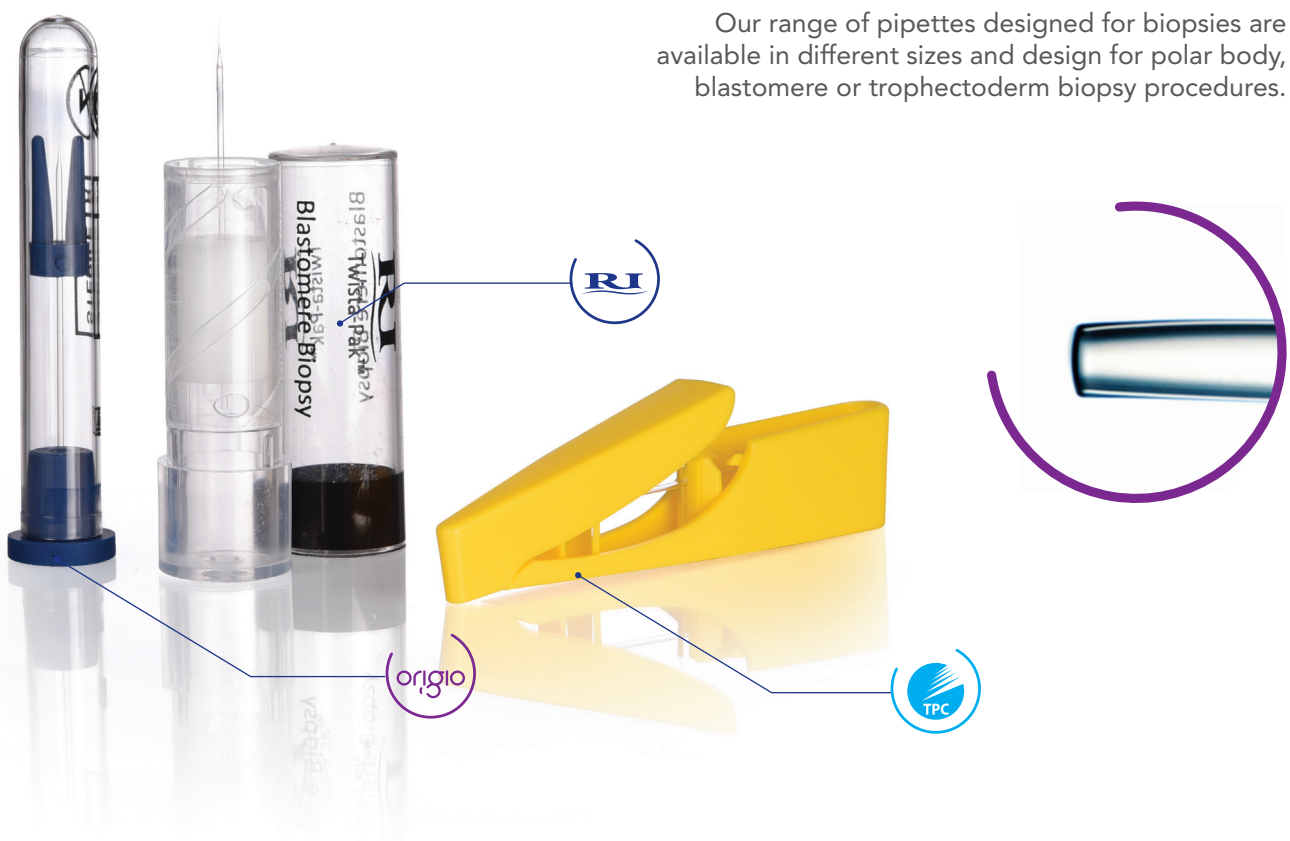

The Stripper® & EZ-Grip®

Tools for Denudation and Handling



Two ranges, tips available in all common sizes

EZ-Grip

- **Precise** - Adjustable aspiration volumes
- **Durable** - Autoclavable, ergonomic design requiring minimal maintenance
- **Versatile** - Fits all plastic denudation and handling tips
- **Fast Identification** - Pipette tip size indicator dial

The Stripper

- **Precise** - Adjustable aspiration volumes
- **Proven** - Tried and tested for decades
- **Compatible** - Fits all plastic denudation and handling tips
- **Varied** - Pipettors available for very small and very large ID tips

EZ-Tip® and Stripper Tips

The close-tolerance plastic tips are flexible, with smooth edges and consistent geometries. Made from unbreakable medical grade plastic, the tips will not scratch Petri dishes.

They are individually packed for sterility and available in multipacks for convenience and economy.



Bulb Pipetting Range

The EZ-Squeeze™ is a medical grade plastic tip fitted with a flexible bulb for a reliable seal and simpler set-up.

EZ-Squeeze tips are available in all common sizes.

EZ-Squeeze

EZ-Squeeze tips are compatible with all EZ-Squeeze handles and Stripper bulb pipettors.

EZ-Squeeze Handle

The autoclavable EZ-Squeeze handle is designed to offer ease of handling and a comfortable working position.



Stripper BP Slimline

The Stripper BP Slimline handle is designed in stainless steel for easy cleaning and autoclaving, and offers simplicity and ease of handling.



Stripper BP

The Stripper BP is specifically designed for increased ergonomics and fits smoothly in your hand for comfort while working.



EZ-Strip®

The EZ-Strip single-use denudation and handling system has been developed for exacting laboratory sterility standards minimising any risk of cross contamination.

It is available in nine sizes, including a long version which is ideal for the handling of gametes in tubes.





Sales & distribution availability of products will vary across markets.

ORIGIO: CE 0086 / Research Instruments: CE 0120 / The Pipette Company: CE 0805

www.origio.com
www.pipetteco.com
www.research-instruments.com



CooperSurgical Companies