



VITRIFICATION BY ORIGIO®

VitriFit™ - Open Carrier

A cool design for optimal storage

- Hexagonal shape for optimal storage efficiency*
- Curved loading area for specimen protection
- Cap and handle designed for perfect fit during storage
- Works for all stages and with all media kits



* Up to 50% improved storage efficiency compared to popular open vitrification carriers

CE
0543

origio


a CooperSurgical Company

Maximize efficiency using minimal storage space with VitriFit


VitriFit is an open system vitrification carrier for storage of oocytes, embryos and blastocysts. It works for all stages and with all media kits.

The hexagonal shape allows VitriFit carriers to be packed tighter during storage than popular vitrification carriers on the market, thereby improving storage efficiency by up to 50%.

The handle and cap are made from the same material and are designed for a perfect fit at all temperatures. The loading area is curved for specimen protection and to ease loading and removal of excess medium to reach minimal volume.

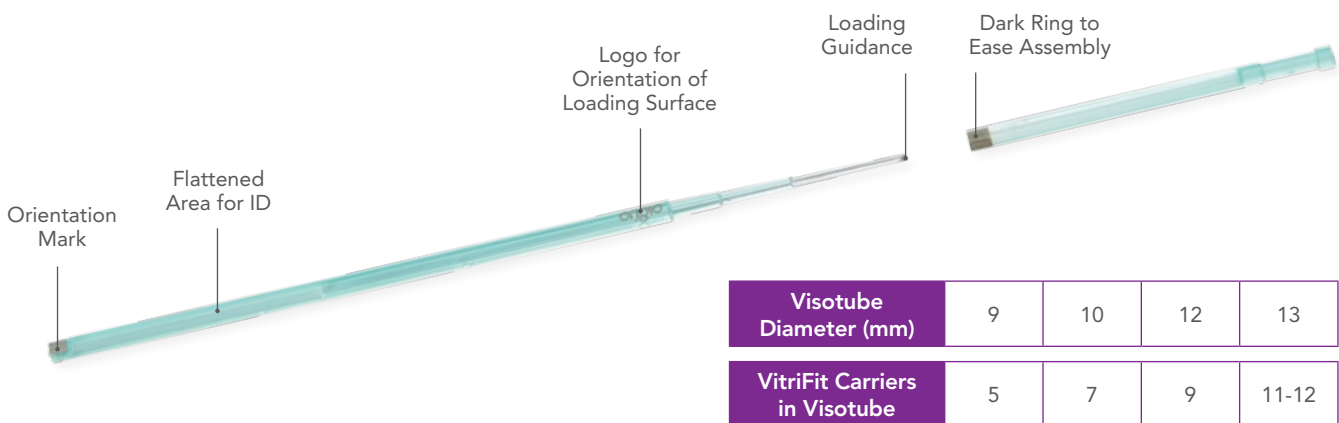

Survival rates for specimen obtained with VitriFit compared to popular open carrier alternatives

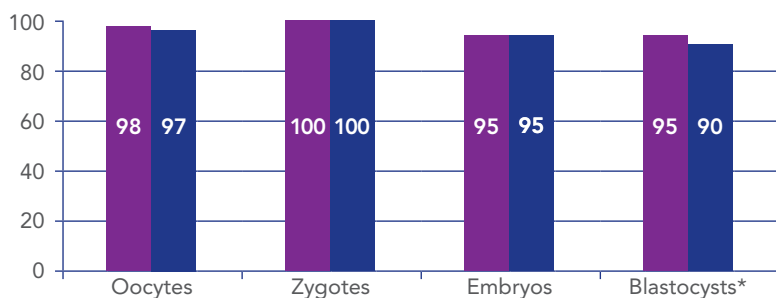

* Mouse tests performed at Embryotech lab inc. Blastocyst survival rate 1 hour post warming

Legend: VitriFit (Purple), Popular Alternative Open Carriers (Blue)

Reference Numbers Box of 20

42782001 - Clear
42792001 - Blue
42802001 - Green
42812001 - Yellow

Learn more about vitrification by ORIGIO

at www.origio.com or contact us at vitrification@origio.com

ORIGIO a/s
 Knardrupvej 2
 2760 Måløv
 Denmark
 Tel: +45 46 79 02 00
www.origio.com



CooperSurgical Companies